




Remember 30-5-15 for injection success



UNLOCK


Hold the pen upright and fully turn the knob to the right to unlock. This starts the process of transferring water to the powder chamber.

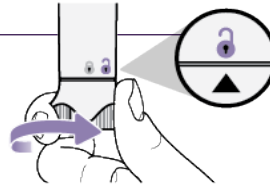

*Turn the knob only once
or the pen may jam.*

SHAKE


Once the gray rod stops moving, shake the pen up and down for **30 SECONDS** to reconstitute.

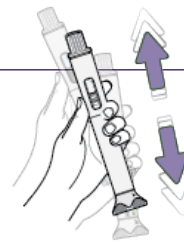

*Ensure the pen is held
upright while shaking to
prevent leaks.*

WAIT


Let the pen sit on a flat surface for **5 MINUTES**. Then, check the inspection window to ensure the powder has completely dissolved.

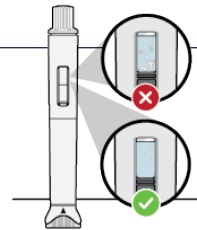

*The powder should be
completely dissolved.
If it isn't, repeat the
previous step.*

TWIST


Hold the pen upright, then twist the cap counterclockwise. After twisting, the cap should easily lift off the pen, and the orange needle guard should have popped up.

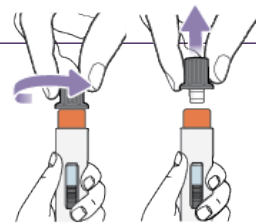

*Hold the pen upright and
gently remove the cap.
Do not force or pull the
cap while twisting.*

INJECT

Identify an injection site on the abdomen 2 inches from the navel or in the upper thigh. Clean the area with an alcohol swab. Inject by placing the pen flat on the skin, gently push down until you hear a click, hold and count for **15 SECONDS**. There is no click to indicate the injection is complete. Lift the pen up and discard safely.

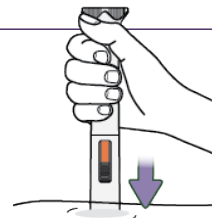

*Do not lift the pen early
or the orange needle
cover will lock.*

Important Safety Information

Indication: NEMLUVIO® (nemolizumab-ilto) is an interleukin-31 receptor antagonist indicated for the treatment of adults with prurigo nodularis. **Contraindication:** Known hypersensitivity to nemolizumab-ilto or to any of the excipients in NEMLUVIO. You are encouraged to report negative side effects of prescription drugs to the FDA. Visit www.fda.gov/medwatch or call 1-800-FDA-1088. **Please see full Prescribing Information, including Patient Information.**

Scan to download
the NEMLUVIO
Instructions for Use



GALDERMA

EST. 1981

©2025 Galderma Laboratories, L.P.
All Rights Reserved. NEMLUVIO is a
registered trademark of Galderma.
US-NLO-2400037 01/25